

For International

**VIVALINE**

 **SHE<sup>N</sup>B**

# VIVALINE

**Comprehensive RF device**

CET·RET·Bipolar·Multi Suction Cup

**Multi Suction RF Technology**

4 multi-suction with 2 channels

**Save your time**

Convenient & effective



# CUPPING SCULPT



## What is VIVALINE?

VIVALINE, as a compound high frequency device, selectively chooses 4 ways of high frequency transfer mode and it includes CET(Capacity Electric Transfer) & RET(Resistive Electric Transfer) & Bi-polar & Multi-suction Cups which bring effective results for skin care, body contouring and pain relief.

## Components



Bipolar  
[1Mhz]



CET + Suction




Multi Suction Cup  
[0.3Mhz]



RET  
[0.5Mhz]

# Multi Suction Effect

## Multi-Cup Principle

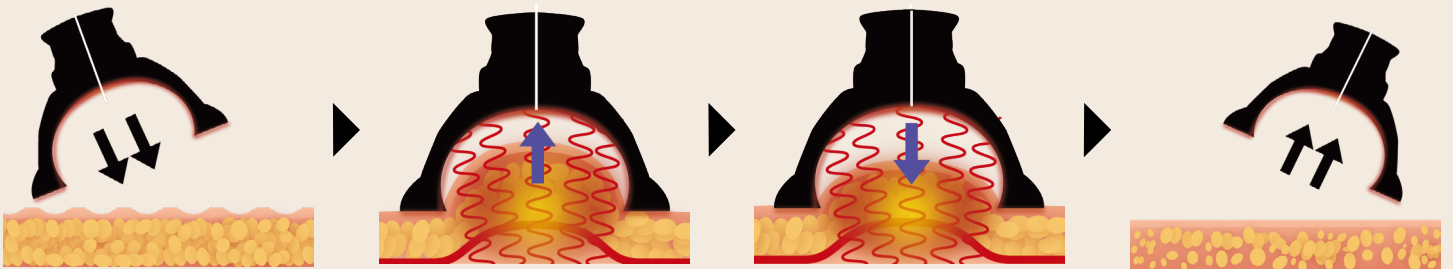

### High Frequency

Skin, body, pain care



### Suction Pulse

Physical therapy using suction & release technology



### Deep Heat

Reach the energy even deeper

## Multi-Cup Advantages

More effective with high frequency and physical therapy

- **Increase blood flow** by negative pressure
- **Increase depth and energy** when transmitting high frequency energy
- **Firming effect** occurred by applying physical energy to the tissue
- **Separation of firm adipose tissue** more effectively
- Accelerates the **regeneration of elastin and collagen** by stimulating the skin tissue
- **Promote the discharge** of lymph glands and decomposed fat cells
- **Relax** the tissue and **Massage** the fascia
- Effective for **cellulite, obesity** and **elasticity** control

# Multi Suction Cup

Cupping Sculpt

---



More effective

---

with high frequency

4 Multi Suction

---

with 2 channels


VIVALINE Patent 10-0922262

An aesthetic medical device that can perform skin care, physical therapy, and obesity management.

In particular, VIVALINE is a multi-function equipped aesthetic medical device, using suction and release by a vacuum pump and radio frequency generation on the contact point of surface of skin simultaneously.

# Indication

## VIVA-Face

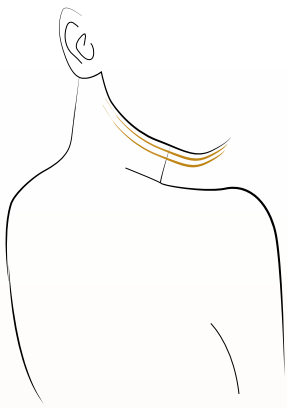

- Lifting · Skin elasticity
- Double Chin · Wrinkle improvement

[Bipolar] + [CET + Suction]


## VIVA-Contouring



- Intra - abdominal fat
- Cellulite · Obesity


[CET + Suction] + [Multi Suction Cup]

## VIVA-Pain

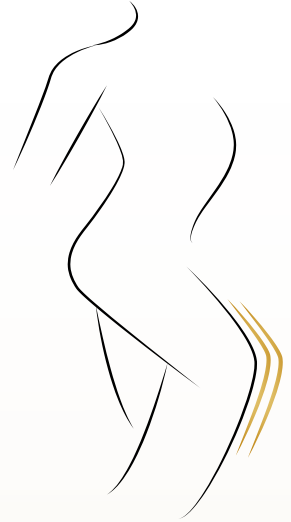

- Arthralgia · Muscle ache
- Sprain · Edema

[Multi Suction Cup] + [RET]

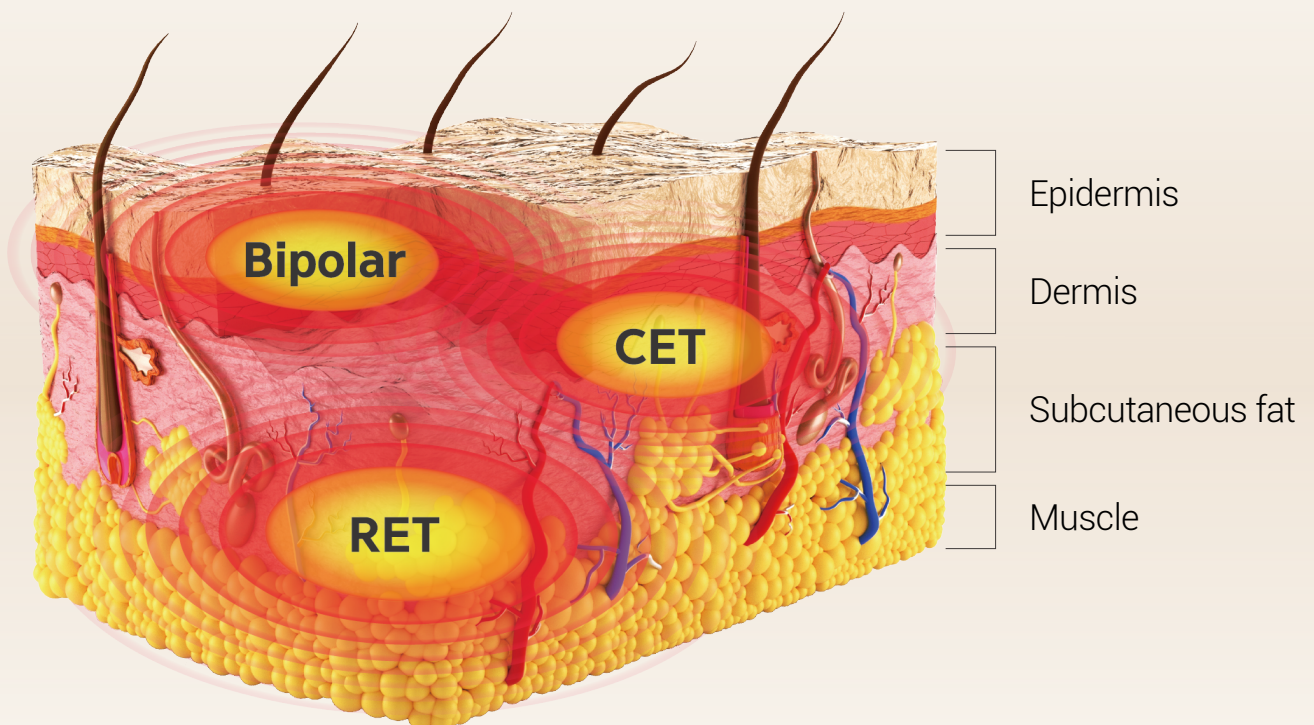

We recorded 4,000 sales world widely together with ex-version of “ V I V A L I N E ” called “ Nemesis ”.





# VIVALINE

<b>Electrical Voltage and Frequency</b>	AC 220V, 50/60Hz
<b>Power consumption</b>	220VA
<b>Weight</b>	4.5kg
<b>Dimension</b>	(W)444.5mm x (L) 400.5mm x (H) 190.1mm
<b>Components</b>	Body, Cable, Handpiece, Suction Cup



[www.shenb.com](http://www.shenb.com)

Tel. +82 2 466 0010 Fax. +82 2 466 5473

SHENB Bldg.,148, Seongsui-Ro, Seongdong-Gu,  
Seoul, Korea